




Complete Tamponade System

ebb is the only dual-balloon tamponade system designed to quickly and effectively treat postpartum hemorrhage.

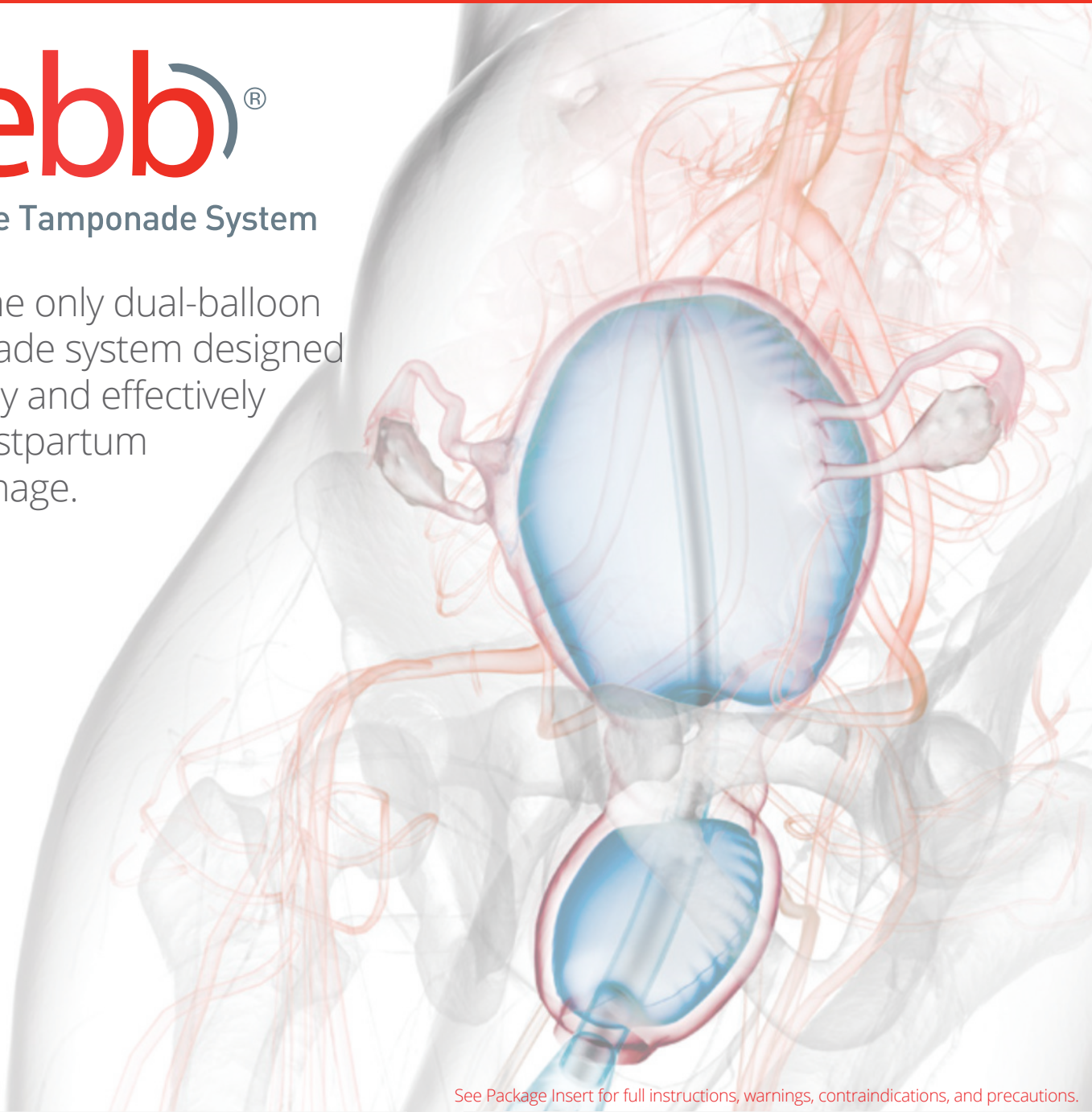

See Package Insert for full instructions, warnings, contraindications, and precautions.

ebb, Designed for Emergency Situations



Conforms to Uterus

The larger (750mL) uterine balloon and malleable material allows ebb to conform to, and treat the uterus.



Rapid Deployment

Simply spike an IV bag and gently squeeze to inflate the balloons in seconds, using the direct-fill ports.



Safety Assurance



Monitor blood loss through the drain and verify lumen patency using the unique irrigation port.




Eliminate Vaginal Packing

The dual-balloon design anchors the device, eliminating the need to pack the vagina and risk of retained pack.


~115,000 mothers with PPH in US each year⁵

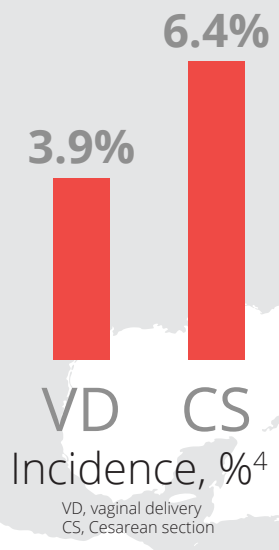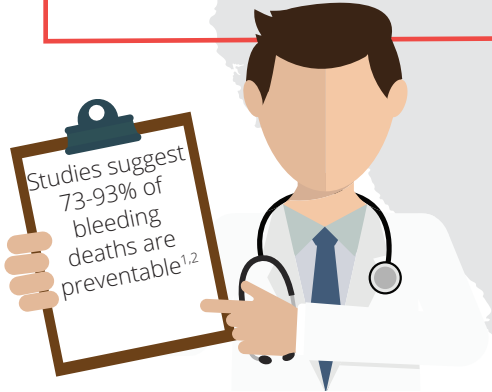

PPH RISK FACTORS^{3,6}

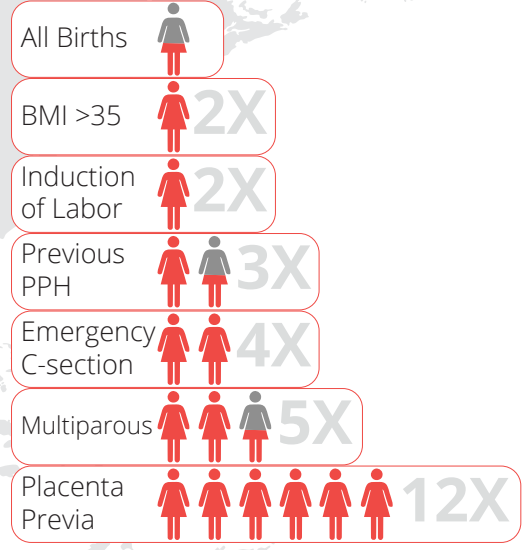

Clinical Evidence of ebb's Efficacy

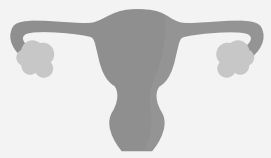

98%
ebb Significantly or Completely Stopped Bleeding

“the [ebb] is very useful in the management of PPH due to uterine atony and abnormal placentation”

Avoided Hysterectomy 92%

45% >500mL


45% patients required >500mL liquid to stop bleeding

90%
Vaginal Balloon Used

Data and quote:
Initial experience with a dual-balloon catheter for the management of postpartum hemorrhage. Dildy GA, Belfort MA, Adair D, et al. Am J Obstet Gynecol. 2014; 210:136.e1-6.

Product Support

ebb is complemented by an offering of training support for tamponade & postpartum hemorrhage.

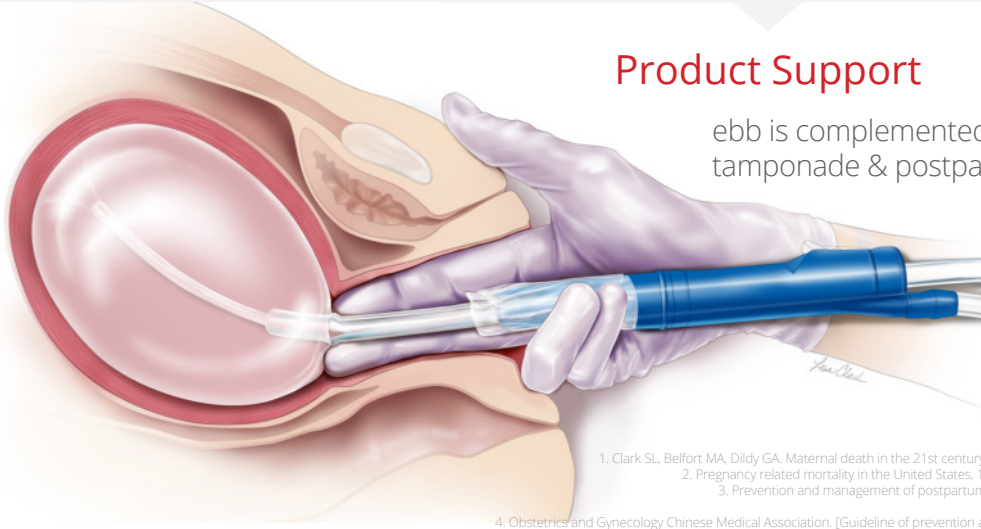

Ordering Information

Product Description	Product Code	Box Quantity
ebb Complete Tamponade System	CTS-1000	1/ea
PPH Trainer (Training Accessory)	052-0083	1/ea

1. Clark SL, Belfort MA, Dildy GA. Maternal death in the 21st century: causes, prevention, and relationship to cesarean delivery. Am J Obstet Gynecol. 2008;199:36e1-36.e5.
 2. Pregnancy related mortality in the United States, 1998 to 2005. Berg CJ, Callaghan WM, Syverson C, Henderson Z. Obstet Gynecol. 2005;106:1228-1234.
 3. Prevention and management of postpartum haemorrhage. Royal College of Obstetricians and Gynaecologists (May 2009 with revisions April 2011). The Journal of Midwifery & Women's Health. 2004; 49(4):283-92.
 4. Obstetrics and Gynecology Chinese Medical Association. [Guideline of prevention and treatment about postpartum hemorrhage (2014)]. Zhonghua Fu Chan Ke Za Zhi 2014; 49: 641-646.
 5. Callaghan WM, Kuklina EV, Berg CJ. Trends in postpartum hemorrhage: United States, 1994 -2006. Am J Obstet Gynecol 2010;202:353.e1-6.
 6. Intrapartum care: care of healthy women and their babies during childbirth; NICE Clinical Guideline (Dec 2014).

Local Distributor:

clinicalinnovations.com
 747 West 4170 South | Murray, Utah USA 84123
 Phone 801.268.8200 | Fax 801.266.7373 |
 Toll Free 888.268.6222

